



Be Smart. Be Sure.

Pick the one that's just right for you.

Only Smartlipo Triplex offers three onboard wavelengths—now including 1440 nm for proficient fat disruption. It's laser lipolysis redefined only from Cynosure.

Smartlipo Triplex™

With 3 wavelengths in a single workstation, you can treat a wider variety of patient concerns more effectively than ever before.

To effortlessly disrupt fatty tissue with less downtime:


Blend **1440 nm** ⊕ **1064 nm**

With 40X greater absorption in fatty tissue – the addition of the 1440 nm wavelength significantly exceeds the performance of all diode lasers (924 nm and 980 nm) on the market.

1064nm creates hemostasis in the blood vessels for less bleeding and bruising

How does the 1440 nm and 1064 nm blend achieve superior results?


- Higher absorption in fatty tissue translates to a more intense photomechanical effect
- Power density concentrated at the tip of the laser fiber creates pressure filled micro-bubbles
- Fat cells are broken up mechanically as well as thermally when micro-bubbles disperse

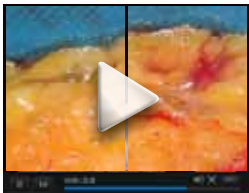

You be the judge.

See for yourself how 1440 nm has the highest absorption of fatty tissue over any other wavelength on the market.
www.cynosure.com/triplex

WAVELENGTH ABSORPTION

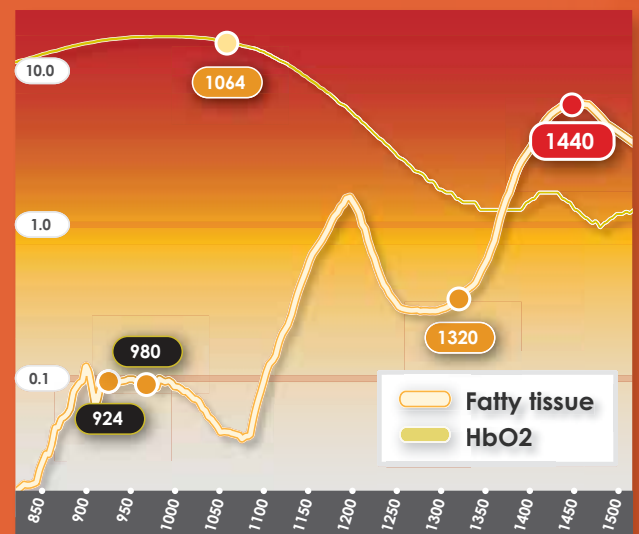

Estimation based on main chromophores in tissue: Human fatty adipose tissue contains approximately 75% lipid and 20% water (Duck, 1990. *Physical Properties of Tissue*, pp 320-320 London: Academic Press)


Smartlipo TriPlex™



To tighten tissue¹ in subdermal layer
with trusted Cynosure wavelengths:

Blend **1320 nm** ⊕ **1064 nm**

CLINICAL OUTCOMES

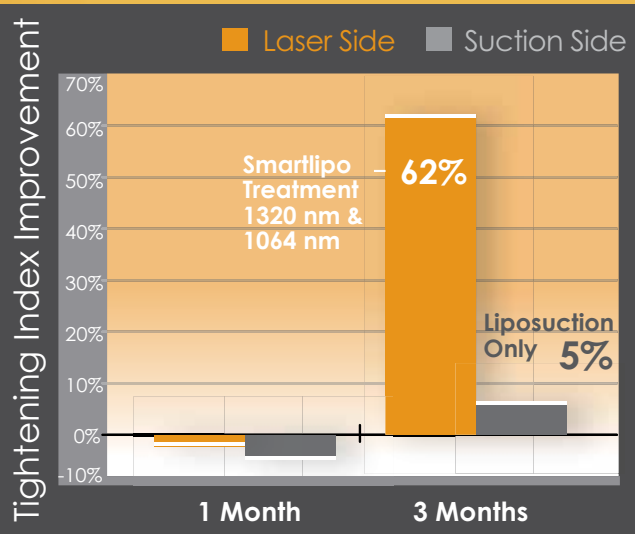

Data courtesy of Barry DiBernardo, MD.
Presented at ASPS 2008.

The proven 1320 nm and 1064 nm blend produces a uniform thermal profile for superior tissue coagulation with tightening that has made Smartlipo so desirable.

Tightening Index number was 57% more than liposuction alone when 1320 nm and 1064 nm are used as reported by Barry DiBernardo, MD.

8 years of clinical studies and experience demonstrate evidence of new collagen production after Smartlipo treatments.

Smartlipo Triplex — laser lipolysis designed just for you — and your patients.

¹ Through tissue coagulation

Choosing the right combination wavelength therapy for optimum body contouring results*

	1320 nm + 1064 nm blend	1440 nm + 1064 nm blend
Personalize your treatment for:		
More Efficient Disruption of Fatty Tissue for faster, deeper treatments		○
Performing high definition treatments		○
Localized heating leaves tissue remaining after aspiration intact for faster recovery		○
Areas of skin laxity where proven tightening results are required	○	
Broad tissue heating for proven skin tightening results through coagulation	○	
Smaller treatment areas with less fat, like the face and neck	○	
Reduced bleeding and bruising in all treatment areas	○	○
Trusted and Reliable Results	○	○

The right wavelength blend and intelligent delivery systems help you deliver the maximum dose, safely.

SmartSense. ThermoGuide. ThermoView. Cynosure's unique, intelligent delivery systems are so smart, they can help you prevent over-treatment, determine your endpoint and minimize complications with integrated temperature and motion sensing.

You'll find all three systems are available for Smartlipo Triplex, the laser lipolysis system that started it all.

Be Smart. Be Sure. Pick the one that's just right for you — Smartlipo Triplex



To find out more about Smartlipo Triplex contact Cynosure today at 800.886.2966 or visit www.cynosure.com

CORPORATE HEADQUARTERS

United States

Cynosure, Inc.
5 Carlisle Road
Westford, MA 01886 USA
Tel: 978.256.4200
Toll-Free: 800.886.2966

* Selection should be based on physician experience or preference

Cynosure is a registered trademark of Cynosure, Inc. • Smartlipo is a trademark of El. En. S.p.A. • Triplex is a trademark of Cynosure, Inc. • LaserBodySculpting is a service mark of Cynosure, Inc. ©2009 Cynosure, Inc. • 921-0177-000 Rev.1 09/2009

